

## INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification : <b>G01N 1/31 // 33/48</b>	<b>A2</b>	(11) International Publication Number: <b>WO 00/36393</b> (43) International Publication Date: 22 June 2000 (22.06.00)
---	-----------	---

(21) International Application Number: PCT/US99/29515


(22) International Filing Date: 14 December 1999 (14.12.99)

(30) Priority Data:  
09/212,367 15 December 1998 (15.12.98) US(71) Applicant: VENTANA MEDICAL SYSTEMS, INC. [US/US];  
3865 North Business Center Drive, Tucson, AZ 85705 (US).(71)(72) Applicants and Inventors: BAUNOCH, David [US/US];  
11038 N. Cloud View, Oro Valley, AZ 85737 (US). JONES,  
Richard [US/US]; 7110 S. Jen Tilly Lane, Tempe, AZ 85283  
(US).(74) Agent: PENN, Amir, N.; McDonnell Boehnen Hulbert &  
Berghoff, 32nd Floor, 300 South Wacker Drive, Chicago,  
IL 60606 (US).(81) Designated States: CA, JP, European patent (AT, BE, CH, CY,  
DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT,  
SE).**Published***Without international search report and to be republished  
upon receipt of that report.*

(54) Title: METHOD AND APPARATUS FOR AUTOMATED REPROCESSING OF TISSUE SAMPLES

## (57) Abstract

A method and apparatus of automatically reprocessing a specimen for microscopic examination is disclosed. Processing of a specimen for microscopic examination involves fixation of the specimen and preparation of the embedded specimen from the fixed specimen. There are instances where, once a specimen has been processed, it is necessary to reprocess the specimen due to contamination of reagents during processing or inadequate fixation. The system automatically reprocesses a specimen by removing residual embedding material from the specimen with a clearing agent, removing the clearing agent with a dehydrating agent, and removing the dehydrating agent with an aqueous fluid.

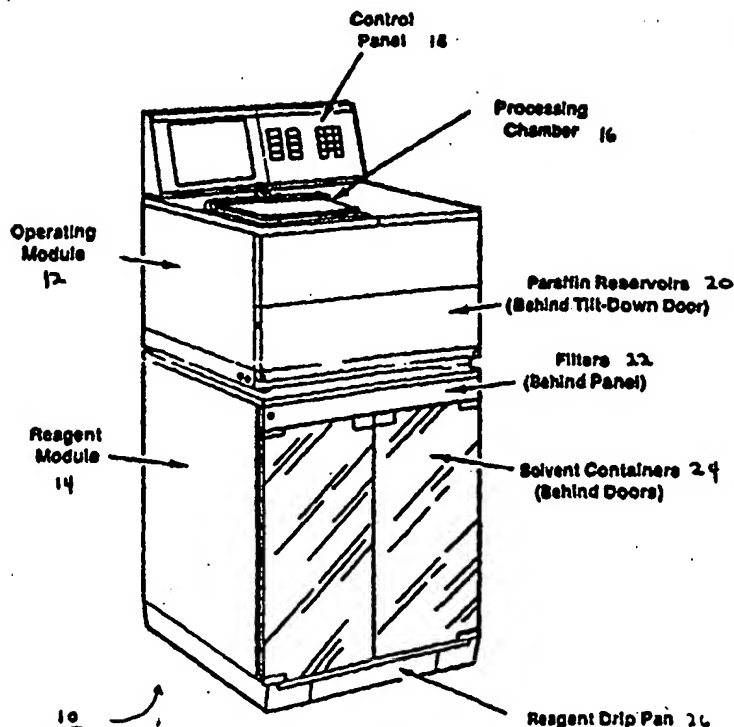

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AL	Albania	ES	Spain	LS	Lesotho	SI	Slovenia
AM	Armenia	FI	Finland	LT	Lithuania	SK	Slovakia
AT	Austria	FR	France	LU	Luxembourg	SN	Senegal
AU	Australia	GA	Gabon	LV	Latvia	SZ	Swaziland
AZ	Azerbaijan	GB	United Kingdom	MC	Monaco	TD	Chad
BA	Bosnia and Herzegovina	GE	Georgia	MD	Republic of Moldova	TG	Togo
BB	Barbados	GH	Ghana	MG	Madagascar	TJ	Tajikistan
BE	Belgium	GN	Guinea	MK	The former Yugoslav Republic of Macedonia	TM	Turkmenistan
BF	Burkina Faso	GR	Greece			TR	Turkey
BG	Bulgaria	HU	Hungary	ML	Mali	TT	Trinidad and Tobago
BJ	Benin	IE	Ireland	MN	Mongolia	UA	Ukraine
BR	Brazil	IL	Israel	MR	Mauritania	UG	Uganda
BY	Belarus	IS	Iceland	MW	Malawi	US	United States of America
CA	Canada	IT	Italy	MX	Mexico	UZ	Uzbekistan
CF	Central African Republic	JP	Japan	NE	Niger	VN	Viet Nam
CG	Congo	KE	Kenya	NL	Netherlands	YU	Yugoslavia
CH	Switzerland	KG	Kyrgyzstan	NO	Norway	ZW	Zimbabwe
CI	Côte d'Ivoire	KP	Democratic People's Republic of Korea	NZ	New Zealand		
CM	Cameroon			PL	Poland		
CN	China	KR	Republic of Korea	PT	Portugal		
CU	Cuba	KZ	Kazakstan	RO	Romania		
CZ	Czech Republic	LC	Saint Lucia	RU	Russian Federation		
DE	Germany	LI	Liechtenstein	SD	Sudan		
DK	Denmark	LK	Sri Lanka	SE	Sweden		
EE	Estonia	LR	Liberia	SG	Singapore		

## METHOD AND APPARATUS FOR AUTOMATED REPROCESSING OF TISSUE SAMPLES

5

### NOTICE REGARDING COPYRIGHT

A portion of the disclosure of this patent document contains matter subject to copyright protection. The copyright owner has no objection to the facsimile reproduction by anyone of the patent disclosure document as it appears in the Patent and Trademark Office files and records but otherwise retains all copyrights whatsoever.

### BACKGROUND OF THE INVENTION

#### A. Field of the Invention

15 This invention relates generally to the fields of histology and cytology, and more particularly relates to a method and apparatus for reprocessing and processing a specimen.

#### B. Description of Related Art

Microscopic examination of specimen samples typically involves examining a slice or a cross-section of the sample. In order to obtain a cross-section, the specimen sample undergoes a process to infiltrate the specimen with a paraffin wax or a wax substitute. Thereafter, the block is embedded and sliced into sections using a microtome.

The method of processing the specimen involves fixation of the specimen and preparation of the infiltrated specimen from the fixed specimen. Fixation of the specimen typically involves immersion, subjecting or exposure of the specimen in a

fixing agent, such as formalin. Preparation of the infiltrated specimen from the fixed specimen is typically a time-consuming, multi-step process requiring dehydration of the fixed specimen with a dehydrant (such as alcohol), clearing of the dehydrant with a suitable clearant (a typical solvent is xylene), and infiltration of the specimen with an infiltrating medium, such as paraffin wax. In addition, the dehydration and clearing steps typically require immersion, subjecting or exposure of the specimen in a graded series of reagents for comparatively long periods of time. The time required for tissue preparation may be on the order of 8 to 12 hours. Examples of tissue preparation are in U.S. Patent No. 3,961,097 entitled "Method of Preparing Tissue for Microscopic Examination" and U.S. Patent No. 4,656,047 entitled "Rapid Method for Cell Block Preparation," both of which are hereby incorporated by reference in their entirety.

Different types of specimens, such as any organelle, cell, cell suspension, tissue section, or tissue specimen, may be infiltrated with a paraffin medium for examination. However, different types of specimens may require different types of procedures to be processed properly. In addition, there may be instances where the specimen may be processed incorrectly, due to contamination of reagents during processing or inadequate fixation. It is typically not until after the specimen has been embedded and sliced that it can be determined whether the specimen has been properly processed. At that point, there are two options: obtain another specimen or reprocess the embedded specimen. If one chooses to reprocess the sample, this involves sequentially immersing, exposing or subjecting the specimen with a series of reagents under controlled conditions. However, this process is very time-consuming and requires a technician to manually proceed through each of the reprocessing steps.

Further, there are instances where a slice or a cross-section of a specimen, after being processed, will be reprocessed for analysis. One instance is ploidy analysis in which tissue sections are cut from the paraffin block, wrapped in a permeable material and reprocessed. The reprocessing steps remove the paraffin using a clearing agent, 5 remove the clearing agent using a dehydrant and remove the dehydrant using an aqueous medium. Nuclei from the specimen are then prepared for DNA analysis using a fluorescent compound.

## SUMMARY OF THE INVENTION

In accordance with a first aspect of the invention, an apparatus for automatically reprocessing a specimen from an infiltrated medium to an aqueous fluid is provided. The apparatus has a processing chamber for holding a specimen, means for regulating flow of fluid to the processing chamber, at least one container of a clearant agent, at least one container of a dehydrant agent and at least one container of an aqueous fluid, the containers of clearant, dehydrant and aqueous fluid being connected to the processing chamber via means for regulating flow of fluid to the processing chamber, and a control device having a processor and a memory device, the processor controlling the means for regulating flow of fluid in order to automatically and sequentially connect the processing chamber with the container of clearant agent, the container of dehydrant agent and the container of aqueous solution in order to reprocess the specimen.

In accordance with a second aspect of the invention, a method for automatically reprocessing a specimen using a specimen reprocessing machine having processor for controlling the exposure of the specimen to a clearing agent, a dehydrating agent and an aqueous fluid is provided. The method includes the step of providing the specimen which is infiltrated with an infiltrating medium, indicating to the specimen reprocessing machine that the specimen is to be reprocessed, exposing the specimen to a clearing agent via the processor to remove the infiltrating medium from the specimen, exposing the specimen to a dehydrating agent via the processor to remove the clearing agent, and exposing the specimen to an aqueous fluid via the processor to remove the dehydrating agent from the specimen.

Accordingly, a goal is to process and reprocess specimens for microscopic

examination. These and other objects, features, and advantages of the present invention are discussed or apparent in the following detailed description.

## BRIEF DESCRIPTION OF THE DRAWINGS

A presently preferred embodiment of the present invention is described herein with reference to the drawings wherein:

FIG. 1 is a front perspective view of the processing and reprocessing system;

5      FIG. 2a is block diagram of the Operating Module of the processing and reprocessing system;

FIG. 2b is block diagram of the Reagent Module of the processing and reprocessing system;

FIG. 3 is a block diagram of the pressure modifier, float valve and processing  
10      chamber in the Operating Module and Reagent Module of Figures 2a and 2b;

FIG. 4 is a front view of the Reagent Module of Figure 1 with the doors removed;

FIG. 5 is a flow chart of the processing of a specimen;

FIG. 6a is a flow chart of the reprocessing of a specimen until introduction of  
15      an aqueous fluid in the specimen and processing of the specimen; and

FIG. 6b is a flow chart of the reprocessing of a specimen, until the step as indicated by the operator, and processing of the specimen.

## DETAILED DESCRIPTION OF PREFERRED AND ALTERNATIVE EMBODIMENTS OF THE INVENTION

20

The processing and reprocessing of tissue is accomplished by sequentially putting the specimen (such as any organelle, cell, cell suspension, tissue section, or tissue specimen) to be processed or reprocessed in contact with, or immersed in, a series of reagents under controlled conditions. The reagents may be divided into three



types: paraffin, solvents and aqueous solvents. The conditions that can be controlled while the tissue is in contact with a reagent can be any combination of heat, pressure, vacuum and agitation.

Referring to Figure 1, there is shown one example of a tissue processing and reprocessing system 10. The tissue processing and reprocessing system may consist of two major components: an Operating Module 12 and a Reagent Module 14. The Operating Module and Reagent Module can be placed side-by-side on a benchtop or stacked for a floor mounted configuration, as shown in Figure 1. The specimen is placed in a processing chamber 16, and reagents are sequentially put into the processing chamber 16 from the solvent containers 24, with excess reagents being collected in the reagent drip pan 26 in case of a malfunction. Paraffin is also introduced into the processing chamber 16 with paraffin reservoirs 20. In this apparatus, the control panel 18 indicates the operation of the system 10 and allows for control of the heat, pressure, vacuum and agitation, which affect the processing chamber 16.

An alternative method of putting the specimen in contact with, or immersed in, the reagents is to have each reagent contained in a separate container, and have a mechanical device, such as a robotic arm and controls, to move the specimen from container to container. With this method, the systems for controlling the heat, pressure, vacuum and agitation can, in any combination, be attached to the individual reagent containers or the device for moving the specimen.

Referring to Figure 2a, there is shown a block diagram of the Operating Module 12. The Operating Module 12 houses the processing chamber 16, the control device 28, Input/Output device 30, and the paraffin oven 32 with three reservoirs 20.

The processing chamber 16 inside the Operating Module 12 connects through tubing to valve 34, such as a rotary valve (or other means for regulating flow of a fluid) to the paraffin reservoirs, and through tubing to the reagents in the Reagent Module. In an alternate embodiment, the means for regulating flow of fluid from the paraffin reservoirs to the processing chamber may be performed by any valve, flap, lid, or plug. The processing chamber has an agitator 36, used when the specimen and reagent require stirring. The agitator 36 may be in the form of a rotating stirring device, a recirculating pump, or any other device that causes the reagent to move with respect to the tissue or the tissue to move with respect to the reagent. In addition, the processing chamber has a pressure sensor 38, used to indicate the pressure in the processing chamber 16 to the pressure modifier 44. As described subsequently, the pressure modifier 44 may be accomplished through mechanical means by applying direct mechanical force to the processing chamber through an aneroid, diaphragm, or other mechanical device. The pressure may also be changed by applying pneumatic pressure or vacuum to the processing chamber (e.g., a compressor (air pump) in the system or an external source of vacuum and/or pressure). This may also be accomplished with a mechanical regulator or by cycling the sources of vacuum or pressure on and off.

The Operating Module 12 also includes the paraffin oven 32. Processing and reprocessing of tissue may include the use of an infiltrating medium such as paraffin. The paraffin is stored in a temperature-controlled container in order to keep the paraffin in a liquid state. The temperature of the paraffin reservoirs 20 can be controlled by applying heat directly to the individual containers or by having the paraffin container(s) in a temperature-controlled chamber (such as an oven 32). The

oven 32 maintains the paraffin in a liquid state so the system can draw the paraffin into the processing chamber 16, allowing it to penetrate the samples. The processing chamber 16 connects to the paraffin reservoirs through a heated rotary valve 34, which facilitates paraffin selection. At the proper time in the processing and reprocessing program, the rotary valve 34 permits paraffin from the selected reservoir to flow into the processing chamber, drawn in under vacuum. During the drain cycle, the valve also selects the proper reservoir for the chamber to empty into.

Referring to Figure 2b, there is shown a block diagram of the Reagent Module 14. Processing and reprocessing may require the use of reagents. The Reagent Module 14 contains reagent containers 40 and is connected to the processing chamber 16 in the Operating Module via a solvent/purge line. In one embodiment, the Reagent Module 14 contains twelve reagent containers: ten solvent containers and two purge containers (as shown in Figure 4). The storage temperature of the solvents typically do not need to be controlled and are therefore stored at room temperature. In this arrangement, there is a means for selecting the specific reagent container 40 to move reagents into the processing chamber 16. In one embodiment, the specific reagent container is selected via a set of two valves, one valve 42 (which is set by the processor 54) in the Reagent Module and the second valve 34 (which is set by the processor 54) in the Operating Module. The valve 42 acts as a means for regulating the flow of fluid (which in a preferred embodiment is a liquid, and in alternate embodiments may include a liquid, gas or both liquid and gas) between the container 40 and the processing chamber. The valve 42 selects which solvent container connects to the fluid line going to the Operating Module. Thereafter, valve 34 selects which of the paraffin lines or solvent/purge line is connected to the processing

chamber. At the proper time in the processing or reprocessing program, the rotary valves 34, 42 (as set by the processor 54) permit only one solvent to flow through the line into the processing chamber, drawn in under vacuum. During the drain cycle, the processor 54 selects the proper setting of valves 34, 42 (as set by the processor 54) for the proper station in order to permit the chamber to empty under pressure. In an alternative embodiment, a valve or other means for regulating the flow of fluid may be mounted on each individual reagent container and connected to a common manifold which connects to the processing chamber. In another embodiment, other means for regulating the flow of liquid between the containers and the processing chamber include any valve, flap, lid, or plug.

The Reagent Module 14 also has another line to the Operating Module 12 that modifies the pressure in the processing chamber 16. The pressure is modified in the tissue processing and reprocessing system via a pressure modifier 44. As shown in more detail in Figure 3, the pressure modifier 44 serves as a means for introducing and extracting reagents from the processing chamber 16. In one embodiment, this is accomplished by using a pump 72 and a series of valves 74, 76, 78, 80. The pump 72 and valves 74, 76, 78, 80 are contained in the pressure modifier 44 to direct and control the pressure and vacuum. A pressure sensor 38 senses pressure or vacuum in the processing chamber 16. The pump 72 in the Reagent Module cycles on and off as needed to lift fluids into the processing chamber 16 and to drain fluids to their containers 40. In a first state of operation, when valves V1, V2 (74, 78) are closed and valves V3, V4 (76, 80) are open, the pump 72 acts to create a vacuum in the processing chamber 16. In this manner, liquids are drawn into the chamber. Air from the processing chamber 16 is sent to the filter 50, which is described subsequently. In

a second state of operation, when valves V3, V4 (76, 80) are closed and valves V1, V2 (74, 78) are open, the pump 72 acts to create a pressure in the processing chamber 16. In this manner, liquids are expelled from the chamber, draining into their respective containers. Air is sent to the processing chamber from a vent. This allows for filling and draining the processing chamber while the specimens remain stationary in the processing chamber. This also permits placing the specimens under vacuum or pressure cycles while immersed in solvents or paraffin to permit thorough infiltration. A pressure differential is created between the storage container and the processing chamber using the pump, to move the reagents. Alternatively, the force of gravity may be used to move the reagents or paraffin. Further, as shown in Figure 3, there is a float valve 82 which prevents reagents from the processing chamber 16 to enter the pressure modifier 44 in the event that the fluid level in the processing chamber 16 is too high.

The Reagent Module 14 has electrical cables for the pump 72, pneumatic valves, 74, 76, 78, 80, rotary valve 42 and blower 48. The Reagent Module 14 features a ventilation system 46 that uses activated charcoal filters to collect solvent fumes before they can escape into the atmosphere. The processing and reprocessing system design reduces the production of fumes. The system handles the fumes from these sources with a built-in ventilating system that filters the air through activated charcoal granules. The system consists of a blower 48 and two filter sections: one for solvent fumes 50, the other for formaldehyde fumes 52. The blower 48 draws air through the Reagent Module 14, up through the filters, and out the back of the unit. This filter system allows operation of the system without the need for a fume hood, external ventilating system, or exhaust fan.

The Operating Module 12 further includes a control device 28. The control device, in one embodiment, may be a general purpose computer. This control device 28 automatically controls and sequences the operation of the heaters, motors, pumps and valves, which are controlled via cables. The control device 28 includes, in a preferred embodiment, a processor 54, and in particular, a Hitachi HD-64180 (Z-80) microcontroller. The control device may also include an electro-mechanical timer, an embedded microprocessor circuit, a programmable logic controller, an external computer, or any combination of the above. The control device 28, in one embodiment, contains memory 56 or other computer readable storage medium, including both random access memory (RAM) 58 and read only memory (ROM) 60 in the form of an erasable programmable read only memory (EPROM). The EPROM contains the system operating program and the text and screen formats for the display. Referring to Appendix A and incorporated herein by reference, there is listed the software having a set of instructions for reprocessing of a specimen. The software is written in Z-80 assembly programming language and is executed on the Hitachi HD-64180 (Z-80) microcontroller.

The control device 28 reads the temperature (via a temperature sensor 39), pressure (via a pressure sensor 38) and the processing chamber fluid level (via a fluid level sensor 37) through the Interface Board 62 and controls the heaters and motors through the Power Board 64. The Power Board 64 contains the drivers 65 for the motors 72, heaters 21, valves 34, 42 and the stirrer 36.

The control device 28 further communicates with the Input/Output device 30 or other user interface. The Input/Output device 30 includes a control panel 18 featuring a monitor 66 such as liquid crystal display (LCD) for displaying menus,

instructions and message. The Input/Output device 30 also includes a keypad 68, such as a numeric keypad and an alpha-numeric keyboard or other means for input such as a mouse. The LCD screen assists in programming and operating the system. Through menus, the screen shows status, guides the operator in writing and running  
5 reprocessing programs, and serves a variety of maintenance functions. During processing, the monitor 66 shows where specimens are in the cycle, the time in each station, the solution in that station, temperature, and vacuum or pressure. The Input/Output device 30 further includes external ports 70 for connections to external devices such as a printer or a phone line.

10 The control device 28, in combination with the Input/Output device 30, gives the system its programming flexibility. The operator can program each of the stations (twelve solvent stations and three paraffin stations) for process time, temperature, vacuum or pressure. The monitor 66 displays all parameters to help the operator while writing the program. A variety of menus give the operator the flexibility of  
15 performing a variety of maintenance and service procedures. A special help function gives on-screen assistance at any time without affecting the present status. The operator can tailor processes to match tissue requirements for different solutions and soaking times as well as a combination of heat, pressure and vacuum.

### **Processing and Reprocessing Cycles**

20 A processing and reprocessing cycle, in one embodiment, consists of filling the processing chamber with a reagent, processing for a programmable amount of time under conditions of controlled temperature, pressure (or vacuum) and agitation. Then draining the reagent back into its storage container. Examples of cycles of the

specimen reprocessing and processing system are the fill cycle, the drain cycle and the process cycle.

### Fill Cycle

As described previously, four pneumatic valves V1, V2, V3, V4 (74, 76, 78, 80) and the pump 72 perform these cycles, all under computer control.

Before the fill cycle, the system checks that the paraffin oven 32, processing chamber 16 and rotary valve block 34 are up to the programmed temperature. The system then vents the processing chamber 16 and calibrates the pressure sensor 38. The solvent rotary valve 42 moves to the proper position for the selected station and the processing chamber rotary valve 34 moves to the closed position for that station. The system then sets the solenoid valves for vacuum and starts the pump 72. This verifies that the processing chamber 16 and pressure modifier 44 do not leak. If the processing chamber 16 cannot maintain vacuum, the solenoid valves cycle five times to clear any contamination from the valve seats. The system makes a second attempt to establish vacuum in the chamber. If the processing chamber 16 still cannot maintain vacuum, the system goes to error standby.

If the system successfully established vacuum, then the system vents the processing chamber 16. The processing chamber rotary valve 34 moves to the open position for the selected station. The system sets the solenoid valves 74, 76, 78, 80 for vacuum and cycles the pump 72 on and off to maintain fill vacuum (4 in. Hg for Stations 1-6, 6 in. Hg for Stations 7-10 and purge stations 15 and 16, 2 in. Hg for paraffin stations 11-13, as shown subsequently in Table 1). The system maintains vacuum until the solution triggers the selected level sensor. The processing chamber rotary valve 34 then closes, and the system vents the processing chamber 16.



Drain Cycle

At the beginning of the drain cycle, the system verifies that the paraffin oven 32, processing chamber 16 and valve block 34 are up to temperature. The system then vents the processing chamber 16 and calibrates the pressure sensor 38 by waiting until  
5 there is no change in pressure for 1/4 second. The system then stores the pressure sensor reading as the ambient pressure.

The solvent rotary valve 42 moves to the selected station if the station is a solvent station, then the processing chamber rotary valve 34 moves to the closed position for that station. The system next sets the solenoid valves for pressure 74, 76,  
10 78, 80 and starts the pump 72. Similar to the fill cycle, it does this to verify that the processing chamber 16 and the pressure modifier 44 do not leak. If the processing chamber 16 cannot maintain pressure, the solenoid valves will cycle five times to clear any contamination from the valve seats. The system then makes a second attempt to establish pressure in the processing chamber. If the processing chamber  
15 still cannot maintain pressure, the system goes to error standby.

The system then releases pressure. The processing chamber rotary valve 34 moves to the open position for the selected station. The system sets the solenoid valves 74, 76, 78, 80 to pressure and the pump 72 starts cycling on and off to maintain drain pressure (1 psi). The system will maintain drain pressure until the processing  
20 chamber 16 can no longer hold pressure, indicating that it is empty (the system senses this by the duration of the pump's running cycle). When the system can no longer hold pressure, it vents the processing chamber 16, then waits five seconds for any remaining fluid to collect in the bottom of the processing chamber 16 and its associated plumbing. The system then sets the solenoid valves 74, 76, 78, 80 to

pressure, and turns the pump 72 on for two seconds to clear the processing chamber 16 and plumbing of any remaining fluid. The system then vents the processing chamber 16 to release any remaining pressure.

### Process Cycle

5 A programming option allows the specimen processing and reprocessing system to alternate pressure and vacuum while processing or reprocessing tissue to enhance the infiltration of the tissue samples. During programming, the operator sets the values: up to seven pounds per square inches of pressure and a vacuum of up to fifteen inches of mercury. Before the system begins the pressure cycle, it vents the  
10 pump while maintaining pressure in the processing chamber. The system does this so that the pneumatic pump starts with no load. After the system starts the pump 72, it sets the solenoid valves 74, 76, 78, 80 to pressure. The cycle runs for 3 minutes at each setting, alternating between vacuum and pressure.

### Processing of Specimen

15 As one example of the specimen processing and reprocessing system, the reagents are arranged in 15 "stations" (3 paraffin stations and 12 solvent stations).

	Station	Reagent	Concentration	Description
	1	Formalin		Fixative
	2	Formalin		Fixative
20	3	Isopropyl Alcohol	70%	Dehydrant
	4	Isopropyl Alcohol	95%	Dehydrant
	5	Isopropyl Alcohol	95%	Dehydrant
	6	Isopropyl Alcohol	100%	Dehydrant
	7	Isopropyl Alcohol	100%	Dehydrant
25	8	Alcohol/Xylene	50/50	Dehydrant
	9	Xylene		Clearant
	10	Xylene		Clearant
	11	Paraffin		
	12	Paraffin		
30	13	Paraffin		
	14	Xylene		Purge Clearant

**Table 1 - Reagent Stations**

Depending on the needs in processing the specimen, any number of stations may be present in the machine. In an alternative embodiment in which the specimen is moved from one container of reagents to the next, there may be 15 such containers, as corresponding to the reagents in Table 1, with as a robotic arm and controls to move the specimen from container to container as necessary.

Referring to Figure 5, there is shown a flow chart of a processing of a specimen. The first step involves fixation of the specimen, as shown at block 84. This typically involves immersing the specimen in Formalin, a fixative. In one embodiment of the invention, the specimen is immersed in, subjected to or exposed to a fixing agent at a station, or a multitude of stations, in a processing machine (*see e.g.*, Stations 1 and 2, as shown in Table 1). However, in processing of the sample, the operator may choose to use both stations, only one station or none of the stations (if the specimen has already been immersed in or exposed to a fixative). The specimen is then dehydrated using a dehydrating agent such as alcohol, as shown at block 86. In one embodiment of the invention, the specimen is dehydrated by immersion in, exposure to or being subjected to a series of alcohol reagents with increasing concentration (*see e.g.*, Stations 3-8, as shown in Table 1). The operator of the processing machine may design a single, or a series, of exposures to alcohol depending on the amount of water contained in the specimen. Thereafter, the specimen is cleared of the dehydrant using a clearing agent, such as xylene (*see e.g.*, Stations 9 and 10, as shown in Table 1), as shown at block 88. Again, depending on processing needs, the specimen may be immersed in, exposed to, or subjected to a

single station or both stations. Thereafter, the specimen is infiltrated with an infiltrating medium such as paraffin (*see e.g.*, Stations 11-13, as shown in Table 1), as shown at block 90.

After the specimen has been processed, the machine should be cleaned in order to minimize contaminants of the reagents upon next use the machine. First, a purge clearant, such as xylene, is used in order to clean the paraffin in the processing chamber and the rotary valve on the processing chamber. Second, a purge dehydrant is used to clean any oily residue, or other contaminants, which may be left in the processing chamber. The purge clearant at station 14 is considered to have more impurities of paraffin and other contaminants than, for example, the clearant at station 9. Further, the purge dehydrant at station 15 is considered to have more impurities of oily residue and other contaminants than, for example, the clearant at station 8.

After the tissue has been processed, it is infiltrated with paraffin, embedded in a paraffin block, and sliced into sections using a microtome. At that point, the operator can determine if the specimen has been processed properly. In one instance, the operator may wish to reprocess the remainder of the sample (*i.e.*, the portion of the specimen which has not been sliced up) until the rehydration of the specimen with an aqueous fluid (to a fixing agent, such as formalin, or to water). If that is the case, the operator indicates, via the control panel 18, that the specimen is to be reprocessed. In addition, if the operator wishes to reprocess a section of the specimen, such as for ploidy analysis, the operator indicates, via the control panel 18, to reprocess the specimen.

#### **Reprocessing of Specimen**

Referring to Figure 6a, there is shown a flow chart of the automatic reprocessing of the tissue until rehydration of the specimen and then processing of the specimen. The system, in one embodiment, may wait until the operator has signaled to reprocess the specimen, as shown at block 92. The software therefore has an initiating routine, waiting until the operator has initiated reprocessing. To reprocess tissue, the infiltrating medium is first removed. Typically, a specimen is not only infiltrated with a medium, but also embedded or encased in the same medium. For example, a specimen may be infiltrated with paraffin, and for ease of slicing, may also be embedded or enveloped with a paraffin shell. To remove the paraffin shell, the operator may simply slice the shell away from the specimen with a knife. Otherwise, the operator may allow the reprocessing machine to remove the shell of paraffin, as shown at block 94. Operator input, via the control panel 18, indicates whether the paraffin station(s) are to be used. This may optionally be done by running processing cycles with one or more paraffin stations. (see e.g., Stations 11-13, as shown in Table 1). The shell of paraffin is removed from the specimen by raising the temperature of the tissue to the melting point of the paraffin that has infiltrated the tissue, as shown at block 96.

If the paraffin stations have already been run, the valve 34 and processing chamber 16 are contaminated with paraffin; therefore, the purge clearant should be used. As described previously, the purge clearant is typically used in cleaning the valve 34 and processing chamber 16 when processing a sample. Thus, the purge clearant is already contaminated with paraffin and may clean the valve 34 and processing chamber 16. As a general matter (even if the paraffin station(s) have not been run), the order of the clearant stations may optionally be from the most

contaminated (with paraffin) to the cleanest. This is due to the fact that in removing the paraffin, the clearing agent may become contaminated. In order to avoid contamination of the "cleaner" clearants, the purge clearant should be used first. Otherwise, the cleaner clearants (such as Stations 9 or 10) would be contaminated with paraffin if used directly after a paraffin step. If that were the case, upon processing of a sample again, the clearant in station 9 or 10 would have to be replaced due to contamination. Therefore, the specimen is immersed in, subjected to or exposed to a purge clearant first, as shown at block 98. (*see e.g.*, Station 14, as shown in Table 1). Typically, the specimen is immersed or exposed to the purge clearant for about 20 minutes with the agitator 36 mixing.

The next step is the removal of residual infiltrating medium from a specimen with "cleaner" clearing agent(s) (an agent that is miscible with the embedding and dehydrating agent), as shown at block 100. This is done by running processing cycles with one or more clearant stations, depending on the needs of reprocessing. (*See e.g.*, Stations 9 and 10, as shown in Table 1). The clearant typically used to remove the paraffin is Xylene.

Again, the valve 34 and the processing chamber 16 may be contaminated with an oily residue left by the clearant. Therefore, the specimen may optionally be immersed in, subjected to or exposed to a purge dehydrant before other dehydrants, as shown at block 102. (*see e.g.*, Station 15, as shown in Table 1). As described previously, the purge dehydrant is typically used in cleaning the valve 34 and processing chamber 16 when processing a sample. Thus, the purge dehydrant is already contaminated. Otherwise, the cleaner dehydrants (such as stations 3-8) would be contaminated with the oily residue if used directly after a clearant step. If that were

the case, upon processing of a sample again, the dehydrants would have to be replaced due to contamination. Therefore, the specimen may be immersed or exposed to the purge dehydrant for about 20 minutes with the agitator 36 mixing.

The next step is the removal of the residual clearing agent by saturating the specimen with a dehydrating agent, as shown at block 104. This step is performed whether or not the purge dehydrant is used. This is accomplished by running processing cycles with dehydrants (typically alcohol) with successively higher concentrations of water in which the specimen is immersed in, subjected to or exposed to dehydrant(s). One or many of the dehydrant stations may be used, depending on the needs of reprocessing. (See e.g., Stations 3-8, as shown in Table 1).

The next step is the removal of the dehydrating agent with an aqueous fluid, as shown at block 106 by the specimen immersed in, exposed to or subjected to the aqueous fluid. The aqueous fluid can be used for storage (such as using an aqueous fluid comprised of water) or used to complete the fixation process prior to the repeating of the specimen processing (such as using an aqueous fluid comprised of a fixative such as formalin).

Optionally, the program may wait to determine if the operator has indicated to process specimen, as shown at block 108. The operator may indicate to process specimen at the beginning of the reprocessing sequence, in the middle of reprocessing, or after reprocessing has completed. Alternatively, the program may immediately processes the specimen without operator input.

The fixed specimen is then processed similar to the process steps of Figure 5. In particular, the refixed specimen is dehydrated using a dehydrant (such as alcohol), as shown at block 110. The dehydrant in the specimen is then replaced using a

clearing agent (such as xylene), as shown at block 112. The clearing agent is then replaced using an infiltrating medium (such as paraffin), as shown at block 114.

Referring to Figure 6b, there is shown a flow chart of an alternate embodiment of the automatic reprocessing of the tissue and then processing of the specimen. The operator, after slicing of the specimen, may be able to determine which step in the previous processing sequence was faulty. For example, if the clearant in the processing sequence was contaminated, upon processing, the clearant may have failed to clear all of the dehydrant, thus leaving the specimen with residual dehydrant. Based on this observation, the operator may enter in the control panel 18 the step to which reprocessing should be done, as shown at block 115. This entry may be stored in a look-up table 59 in RAM 58, so that upon reprocessing, the software may determine which step to reprocess to. Alternatively, the entry in the look-up table may be the step in the processing sequence which was faulty. In the example given above, the step would be the clearing step. In this manner, the reprocessing program may read the entry in the look-up table 59, and stop the reprocessing either at the faulty step or the step prior to the faulty step.

Similar to Figure 6a, the system waits until the operator has signaled to reprocess the specimen, as shown at block 116. The operator may allow the reprocessing machine to remove the shell of paraffin, as shown at block 118. The shell of paraffin is removed from the specimen by raising the temperature of the tissue to the melting point of the paraffin that has infiltrated the tissue, as shown at block 120. If the paraffin stations have already been run, the order of the clearant stations is from the most contaminated (with paraffin) to the cleanest. Therefore, the purge



clearant is used first, as shown at block 122. (*see e.g.*, Station 14, as shown in Table 1).

The next step is the removal of residual infiltrating medium from the specimen with a clearing agent (an agent that is miscible with the embedding and dehydrating agent), as shown at block 124. Thereafter, the look-up table 59 is examined to determine whether the clearing step is the last or final step in the reprocess, as shown at block 126. If so, the program then determines if the operator has indicated to process the specimen, as shown at block 128. If so, the program goes to the infiltrating step, as shown at block 146. If the clearing step is not the last or final step in the reprocess, the purge dehydrant is used before other dehydrants, as shown at block 130. (*see e.g.*, Station 15, as shown in Table 1).

The next step is the removal of the residual clearing agent by saturating or exposing the specimen with a dehydrating agent, as shown at block 132. This step is performed whether or not the purge dehydrant is used. One or many of the dehydrant stations may be used, depending on the needs of reprocessing. (*See e.g.*, Stations 3-8, as shown in Table 1).

Thereafter, the look-up table 59 is examined to determine whether the dehydrating step is the last or final step in the reprocess, as shown at block 134. If so, the program then determines if the operator has indicated to process the specimen, as shown at block 136. If so, the program goes to the clearing step, as shown at block 144. If the dehydrating step is not the last step in the reprocess, the next step is the removal of the dehydrating agent with an aqueous fluid, as shown at block 138. Optionally, the program then waits to determine if the operator has indicated to process specimen, as shown at block 140; otherwise, processing begins without

operator input. The operator may indicate to process specimen at the beginning of the reprocessing sequence, in the middle of reprocessing, or after reprocessing has completed. The fixed specimen is then processed similar to the process steps of Figure 5. In particular, the refixed specimen is dehydrated using a dehydrant (such as alcohol), as shown at block 142. The dehydrant in the specimen is then replaced using a clearing agent (such as xylene), as shown at block 144. The clearing agent is then replaced using an infiltrating medium (such as paraffin), as shown at block 146.

From the foregoing detailed description, it will be appreciated that numerous changes and modifications can be made to the hardware and software aspects of the invention without departure from the true spirit and scope of the invention. For example, the present invention is not dependent on any specific type of computer architecture or type of protocol. This true spirit and scope of the invention is defined by the appended claims, to be interpreted in light of the foregoing specification.

## Appendix A

```

*****
*
*      REVERSE RUN
*
5 *****

REVRUN:      CALL  DISREV

10 REVR2:CALL  GETKEY
      CP      KBF1      ;START
      JP      Z,RUNREV
      CP      KBF2      ;EDIT
      JP      Z,EDITREV
15      CP      KBF3      ;REPROCESSING PROGRAM
      JP      Z,REPCHG
      CP      KBF4      ;CANCEL
      JP      Z,MAINTN
      JR      REVR2

20      RET

*****
*
*      DISPLAY PROGRAM - REVERSE
*
25 *****

DISREV:      LD      A,(REVSTEP)
30      PUSH  AF
      LD      A,(STANUB)
      PUSH  AF

      LD      A,81
35      LD      (SCRNUB),A
      CALL  DISPSR

      LD      A,2
      LD      (DISROW),A
40      LD      A,3
      LD      (DISCOL),A
      CALL  PTSET      ;POINT AT PROGRAM

      XOR     A
45      DISRE1: INC     A
      PUSH  AF

      LD      (REVSTEP),A
      CALL  SOREV
50      CALL  GETRRS
      CALL  BLKDAT      ;BLANK DSPDAT
      LD      HL,DSPDAT ;BUFFER AREA
      LD      (HL),00
      INC     HL
55      LD      A,(STANUB)
      CALL  BINASC

```

```

INC    HL
CALL   LINEDT           ;LINE AT A TIME

5      LD    A,(DISROW)
INC    A
LD     (DISROW),A
LD     A,2
LD     (DISCOL),A
CALL   PTSET           ;POINT AT PROGRAM
10     LD    HL,DSPDAT
CALL   MESOUT

15     POP    AF
CP     15
JR     NZ,DISRE1

20     XOR    A
LD     (EDVAR),A       ;ENABLE PROCESS TIME ROUTINE
CALL   REVTP           ;DISPLAY TOTAL PROCESS TIME
CALL   DISS           ;DISPLAY STIR SPEED
CALL   DISL           ;DISPLAY LEVEL SENSOR
CALL   DSPRPN         ;DISPLAY REPROCESSING NAME

25     POP    AF
LD     (STANUB),A
POP    AF
LD     (REVSTEP),A
RET

30     *****
*
*     REVERSE RUN - RUN
*
35     *****

RUNREV: XOR    A
LD     (BEEPFLG),A     ;TURN BEEPER OFF

40     IN0    A,(PORT1C)
RES    5,A             ;DISABLE ALARM RELAY
OUT0   (PORT1C),A

45     CALL   ALLOFF
CALL   ERRSAV

LD     A,(FLAG3)       ;CHECK IF POT FULL
BIT    6,A
JP     NZ,CUSERV       ;IF SO GO TO USER SERVICE

50     LD     A,84
LD     (SCRNUB),A      ;SCREEN 1.2.4
CALL   FUNCT           ;POWER FAIL REENTRY

55     LD     A,1
LD     (REVSTEP),A

```

	CALL	ERRCLR	;CLEAR ALL ERRORS.
	LD	HL,FPCDAT	;SET HOLD TIME TO CURRENT TIME
	LD	DE,HDATE	
5	LD	A,4	
	CALL	COPYN	
	LD	HL,FPCTIM	
	LD	DE,HTIME	
10	LD	A,4	
	CALL	COPYN	
	CALL	ACTREV	;GET FIRST ACTIVE STATION.
	OR	A	;IS ANY ACTIVE?
15	JP	NZ,RRUN2	;EXIT IF NO STATIONS ACTIVE
	LD	A,(REPNUB)	
	OR	A	
	JP	NZ,RRUN14	
20	JP	MAINTN	
	RRUN2:CALL	ACTREV	;GET NEXT ACTIVE STATION.
	LD	A,(REVSTEP)	
25	LD	(LASTST),A	
	CALL	STTITL	;INITIALIZE STATS.
	CALL	ITLVRP	;INIT LEVEL STAT
	CALL	SOREV	
30	LD	A,(STANUB)	
	CALL	GETRRS	
	CALL	SETSTA	
	CALL	FINISH	
35	XOR	A	
	LD	(DSPRMS),A	
	CALL	SPTIME	
40	RRUN3:LD	A,1	;RUN STATE FILL
	LD	(BSTATE),A	
	CALL	DISREV	
45	LD	A,(REVSTEP)	
	CALL	PGMPTR	
	LD	A,(FLAG3)	;THIS IS FOR REENTRY
50	BIT	6,A	
	JR	NZ,RRUN6	;IF POT FULL DO NOT RESTART STATION
	LD	A,(FLAG1)	
	SET	0,A	;INFILTRATION FLAG
55	LD	(FLAG1),A	
	CALL	LIDCK	

	CALL	GETRRS	;GET STATION DATA
	CALL	SETSTA	;SET STATION PARAMETERS.
	CALL	FILPOT	;FILL TISSUE POT
5			
	LD	A,(STIRSP)	
	LD	(STIRS1),A	
	CALL	SPEED	
10			
	LD	A,(STANUB)	;CLEAR PURGE FLAG IF PURGE STATIONS RUN
	CP	14	
	JR	NZ,RRUN4	
	LD	A,(PRGFLG)	
	RES	0,A	
15			
	LD	(PRGFLG),A	
	JR	RRUN5	
	RRUN4 CP	15	
	JR	NZ,RRUN5	
20			
	LD	A,(PRGFLG)	
	RES	1,A	
	LD	(PRGFLG),A	
	RRUN5 LD	A,6	;INIT V/P FLAGS.
25			
	LD	(VPTIM),A	
	CALL	STATIT	;INITIALIZE POT STAT.
	RRUN6:CALL	LIDCK	
30			
	LD	A,2	;RUN STATE PROCESSING
	LD	(BSTATE),A	
	CALL	DISREV	
35			
	XOR	A	
	LD	(DSPRMS),A	;BLANK RUN MESSAGE.
	CALL	ERRSAV	
40			
	XOR	A	
	LD	(ESTATE),A	;RESET ERROR STATE
	LD	A,(STANUB)	
45			
	CALL	GETRRS	
	CALL	SETSTA	
	LD	A,(REVSTEP)	
	CALL	PGMPTR	
50			
	CALL	CLRRX	;CLEAR ALL KEYS ENTERED
	CALL	STOPST	;ENABLE STOP KEY
	LD	A,(FLAG4)	
55			
	RES	3,A	;ENABLE VALVES
	LD	(FLAG4),A	

```
LD      A,(DSPFLG)
RES     0,A          ;DISPLAY ALL FIRST TIME
LD      (DSPFLG),A
```

5 \*\*\*\*\*

```
RRUN7:CALL INPUT          ;TEST FOR STOP

CALL STIRON
10 CALL POTON
CALL VPRUN

CALL DSPRUN
CALL STATCHK
15 CALL LIDCK

LD      A,(CLKFLG)
BIT     0,A          ;TIMED OUT YET?
JP      Z,RRUN7
20
```

\*\*\*\*\*

```
CALL ALLOFF
CALL STOPUS          ;DISABLE STOP KEY
25

LD      A,(STANUB)
CALL ERRSAV
CALL SVSTAT

30 LD      A,(REVSTEP)
INC     A
LD      (REVSTEP),A

CALL ACTREV          ;GET NEXT ACTIVE STATION.
35 OR     A
JR      Z,RRUN10

RRUN8:LD      A,(LASTST)  ;RECALL LAST
40 LD      (REVSTEP),A

CALL DISREV

CALL SOREV
LD      A,(STANUB)
45 CALL GETRRS

LD      A,(REVSTEP)
CALL PGMPTR          ;POINT TO STATION.

50 LD      A,3
LD      (BSTATE),A    ;RUN STATE DRAIN

CALL DRAIN

55 LD      A,1          ;MSG1 LOC.
LD      (DSPRMS),A
CALL RUNMSG          ;DISPLAY DRAIN MESS
```

```

RRUN9 CALL INPUT                                ;CHECK FOR STOP

5      LD  A,(CLKFLG)
      BIT  1,A
      JR   Z,RRUN9

      CALL ERRSAV

10     LD  A,(REVSTEP)
      INC  A
      LD   (REVSTEP),A

      JP   RRUN2

15     RRUN10: CALL ALLOFF                        ;PROCESS COMPLETE

      CALL DISREV

20     LD  A,(REPNUB)
      OR   A
      JR   NZ,RRUN13

      LD  A,(LASTST)                            ;RECALL LAST
      LD  (REVSTEP),A
      CALL SOREV
      LD  A,(STANUB)
      CALL GETRRS
      CALL SETSTA

30     LD  A,8                                  ;SCREEN 1.2.4
      LD  (SCRNUB),A                            ;POWER FAIL REENTRY
      CALL FUNCT
      LD  A,(REVSTEP)
      CALL PGMPTR

35     CALL STIRON

      LD  A,5
      LD  (BSTATE),A                            ;RUN STATE PROCESS COMPLETE

      LD  A,(REVSTEP)
      CALL PGMPTR                            ;POINT TO STATION.

45     LD  A,6                                  ;PROCESSING COMPLETE.
      LD  (DSRMS),A
      CALL RUNMSG

      RRUN11 CALL GETKEY
50     CP   KBF1                                ;DRAIN
      JR   Z,RRUN12
      JR   RRUN11

      RRUN12: CALL ALLOFF

55     LD  A,(STANUB)
      CALL SVSTAT

```



```

CALL ERRSAV

CALL DRAIN

5      CALL ERRSAV      ;SAVE ERRORS
      XOR  A
      LD   (BSTATE),A   ;DEACTIVATE RUN

      LD   A,(PRGFLG)    ;CHECK IF PARAFFINS
10     AND  03H          ;HAVE RUN.
      LD   (PRGFLG),A

      LD   A,(FLAG2)
      RES  7,A
15     LD   (FLAG2),A    ;CLEAR LID OPEN FLAG

      JP   MAINTN

RRUN13: LD   A,(LASTST) ;DRAIN AND JUMP TO PROCESSING PROGRAM
20      LD   (REVSTEP),A ;RECALL LAST

      CALL DISREV

      CALL SOREV
25     LD   A,(STANUB)
      CALL GETRRS

      LD   A,(REVSTEP)
      CALL PGMPTR        ;POINT TO STATION.
30     LD   A,6
      LD   (BSTATE),A    ;RUN STATE DRAIN

      CALL DRAIN

35     RRUN14: CALL ERRSAV ;SAVE ERRORS
          XOR  A
          LD   (BSTATE),A ;DEACTIVATE RUN

          LD   A,(PRGFLG) ;CHECK IF PARAFFINS
40         AND  03H        ;HAVE RUN.
          LD   (PRGFLG),A

          LD   A,(FLAG2)
          RES  7,A
45         LD   (FLAG2),A  ;CLEAR LID OPEN FLAG

          LD   A,(REPNUB)
          LD   (PGMNUB),A
50         CALL GETPGM

          CALL FINISH      ;COMPUTE IMEDIATE FINISH TIME
          XOR  A
          LD   (STANUB),A

55        JP   RUN1

```

\*\*\*\*\*

\*  
\* GETS NEXT ACTIVE STATION AND PUTS IT IN REVSTEP. 00H IF NOT. - REV  
\*

\*\*\*\*\*

5  
ACTREV: LD A,(REVSTEP)  
PUSH AF  
  
10 ACTR1: CALL SOREV  
LD A,(STANUB)  
  
CALL STONCK  
OR A ;IS STATION OFF?  
JR NZ,ACTR3  
  
15 LD A,(REVSTEP)  
CP 15 ;IS NUMBER BEYOND 15?  
JP P,ACTR2  
INC A ;ADVANCE STATION.  
20 LD (REVSTEP),A  
JR ACTR1  
  
ACTR2: POP AF  
LD (REVSTEP),A  
25 XOR A  
RET

ACTR3: POP AF  
RET

30

\*\*\*\*\*

\*  
\* OUTPUT TOTAL PROCESS TIME - REV  
\*

35 \*\*\*\*\*

REVTPT: PUSH AF  
PUSH DE  
PUSH HL  
  
40 LD A,(EDVAR)  
OR A  
JR NZ,REVT2  
  
45 LD HL,FDATE  
LD DE,SDATE  
LD A,6  
CALL COPYN  
  
50 LD DE,CTIM1 ;ZERO ACCUM  
LD HL,ZERO  
CALL KOPY  
  
XOR A  
55 LD (STAON1),A ;INIT STATION BITS.  
LD (STAON2),A

```
LD      A,(STANUB)
CALL    SAVRRS      ;SAVE STADAT
PUSH    AF

5      LD      A,15
      LD      (STANUB),A

      REVT1: CALL  GETRRS      ;GET NEXT STATION
      LD      BC,TIME1
10     LD      HL,TMPR1
      CALL    CFLPT      ;CONVERT TO FLOATING

      CALL    STAON      ;CHECK IF STATION ON.
      LD      IY,CTIM1
15     LD      BC,TMPR1
      LD      HL,CTIM1
      CALL    TIMADD      ;ACCUM TIME

      LD      A,(STANUB)
20     DEC     A
      LD      (STANUB),A
      JR      NZ,REVT1      ;TEST IF 15 STATION

      LD      HL,PTIME
25     LD      BC,CTIM1      ;SAVE TOTAL PROCESS TIME
      LD      A,20H
      CALL    CASC1

      LD      A,18      ;SET POINTER TO STATION
30     LD      (DISROW),A
      LD      A,10
      LD      (DISCOL),A
      CALL    PTSET
      CALL    BLKDAT

35     LD      HL,DSPDAT
      LD      DE,PTIME      ;OUTPUT PROCESS TIME
      CALL    SFTIME
      LD      (HL),1AH

40     LD      HL,DSPDAT
      CALL    MESOUT

      POP     AF
45     LD      (STANUB),A
      CALL    GETRRS      ;RESTORE STADAT

      REVT2 POP     HL
      POP     DE
50     POP     AF
      RET
```

```
*****
*
55  *      REVERSE RUN - EDIT
*
*****
```

```

EDITREV: LD  A,82
          LD  (SCRNUB),A
          CALL FUNCT          ;DISPLAY SCREEN
5
          LD  A,(FLAG1)
          RES 4,A
          LD  (FLAG1),A      ;CLEAR EDITOR RANGE ERROR

10        LD  A,1
          LD  (REVSTEP),A
          XOR  A
          LD  (EDVAR),A

15  RRE1: LD  A,(REVSTEP)
          CALL SOREV
          CALL GETRRS

          LD  A,(FLAG1)
20        SET 7,A
          LD  (FLAG1),A      ;SET REV VIDEO FLAG

          LD  A,(REVSTEP)
          CALL EDLINE

25  RRE2: CALL GETKEY
          CP   KBUP           ;REV FIELD
          JP   Z,RRFREV
          CP   KBDN           ;ADV FIELD
          JP   Z,RRFFWD
30        CP   KBF1           ;F1 DONE
          JP   Z,REVRUN
          CP   KBF2           ;F2 SOLUTION EDIT
          JP   Z,REVSOL
          CP   KBF3           ;F3 CHANGE STIR SPEED
          JP   Z,REVSPD
          CP   KBF4           ;F4 CHANGE LEVEL SENSOR
          JP   Z,REVLEV
          CP   0DH             ;ENTER
40        JP   Z,RRDAT
          CP   KBCL           ;CLEAR ON/OFF
          JR   Z,RONOFF
          AND  07FH
          SUB  030H
          CP   0AH            ;NUMBER
45        JP   P,RRE2
          JP   REVNUB

RONOFF   CALL ONOFF
50        CALL GETBAC
          CALL SAVRRS
          LD  A,(REVSTEP)
          CALL EDLINE
          JR   RRE2
55

```

\*\*\*\*\*

\*

```

*      CHANGE STIR SPEED
*
*****

5  REVSPD:  LD    A,(STIRSP)
          INC    A
          CP     11
          JP     M,REVSPD1
          LD     A,0

10  REVSPD1  LD     (STIRSP),A
          CALL   DISSS
          JP     RRE2

15  *****
*
*      CHANGE LEVEL SENSOR
*
*****

20  REVLEV:  LD     A,(LEVFLG)
          XOR    002H          ;TOGGLE
          LD     (LEVFLG),A
          CALL   DISLS
25  JP      RRE2

*****
*
*      CHANGE REPROCESSING PROGRAM
*
30  *****

REPCHG:  LD     A,(REPNUB)
          INC    A
35  CP     12
          JP     M,REPCHG1
          LD     A,0

40  REPCHG1  LD     (REPNUB),A
          CALL   DSPRPN
          JP     REVR2

*****
*
45  *      INPUT NUMBER INTO VARIABLE POINTED
*      TO BY REVSTEP AND EDVAR
*
*****

50  REVNUB:  PUSH  AF
          CALL   PUTNUB

          LD     A.20H
          CALL   FILVAR          ;FILL VARIABLE WITH BLANKS

55  POP     AF
          ADD    A.30H          ;CONVERT TO ASCII.

```

```

PUSH HL
POP IY
LD (IY+5),A ;LOAD FIRST CHAR.

5 LD A,C
CP 5
CALL Z,EDDP ;CHECK IF TIME VARIABLE.
CP 1
CALL Z,EDDP ;CHECK IF TIME VARIABLE.

10 LD A,(REVSTEP)
CALL EDLINE ;OUTPUT VARIABLE.

REDN1:CALL GETKEY
15 CP KBF1 ;CANCEL
JP Z,RRE1
CP 0DH ;ENTER
JP Z,RRDAT
CP KBCL ;CLEAR
20 JR Z,AONOFF
AND 07FH
SUB 030H
CP 0AH ;NUMBER
JP P,REDN1
25 DEC B
JR Z,REDN2
ADD A,30H
CALL SHFTIN ;SHIFT IN NEW NUMBER.
LD A,C
30 CP 5
CALL Z,EDDP ;CHECK IF TIME VARIABLE.
CP 1
CALL Z,EDDP ;CHECK IF TIME VARIABLE.

35 LD A,(REVSTEP)
CALL EDLINE ;OUTPUT VARIABLE.
JR REDN1

REDN2:INC B
40 JR REDN1

AONOFF CALL ONOFF
CALL GETBAC
CALL SAVRRS
45 LD A,(REVSTEP)
CALL EDLINE
JR REDN1

*****
50 *
* SAVES THE FIELD
*
*****

55 RRDAT:CALL RNGCHK
CP 0AAH ;TEST IF IN RANGE.
JP NZ,RRD1

```

```

      CALL GETRRS
      LD  A,(FLAG1)
      SET 4,A
      LD  (FLAG1),A      ;SET EDITOR RANGE ERROR
5      JP  RRE1

RRD1  CALL SAVRRS
      LD  A,(REVSTEP)
      CALL OUTLIN
10     LD  A,(FLAG1)
      RES 4,A
      LD  (FLAG1),A      ;CLEAR EDITOR RANGE ERROR
      JP  RRFFWD

15     *****
      *
      *   REV FIELD
      *
      *****

20     RRFREV: LD  A,(FLAG1)
      RES 7,A
      LD  (FLAG1),A      ;CLEAR REV VIDEO FLAG

25     BIT 4,A      ;RETURN IF EDITOR RANGE ERROR
      JP  NZ,RRE1

      LD  A,(REVSTEP)
      CALL EDLINE
30     CALL REVTPY      ;DISPLAY TOTAL PROGRAM TIME

      LD  A,(EDVAR)
      OR  A      ;FIRST VARIABLE IN STA?
      JR  Z,RRF1
35     DEC  A
      LD  (EDVAR),A
      JP  RRE1

RRF1:  CALL SAVRRS
40     LD  A,(REVSTEP)
      CALL OUTLIN
      CP  1      ;FIRST STATION?
      JR  Z,RRF2

45     DEC  A      ;ACCESS PREVIOUS LINE
      LD  (REVSTEP),A
RRF4:  CALL OFFRCK
      OR  A
      JR  NZ,RRF3
50     XOR  A
      LD  (EDVAR),A
      JP  RRE1

RRF2:  LD  A,15
55     LD  (REVSTEP),A
      JR  RRF4

```

```

RFF3: LD  A,3                ;RAP AROUND
      LD  (EDVAR),A
      JP  RRE1

5  *****
  *
  *   ADVANCE FIELD
  *
  *****

10  RRFFWD: LD  A,(FLAG1)
      RES  7,A
      LD  (FLAG1),A          ;CLEAR REV VIDEO FLAG

15  BIT  4,A                ;RETURN IF EDITOR RANGE ERROR
      JP  NZ,RRE1

      LD  A,(REVSTEP)
      CALL EDLINE
20  CALL REVTP             ;DISPLAY TOTAL PROGRAM TIME

      RAF0 LD  A,(EDVAR)
          CP  3                ;FIRST VARIABLE IN STA?
          JR  Z,RAF1

25  CALL OFFRCK
      OR  A
      JR  Z,RAF1

      LD  A,(EDVAR)
30  INC  A
      LD  (EDVAR),A
      JP  RRE1

      RAF1: CALL SAVRRS
          LD  A,(REVSTEP)
          CALL OUTLIN
          CP  15                ;FIRST STATION?
          JR  Z,RAF2
          INC A                ;ACCESS PREVIOUS LINE
40  LD  (REVSTEP),A
      JR  RAF3

      RAF2: LD  A,1
          LD  (REVSTEP),A

45  RAF3: XOR  A                ;RAP AROUND
          LD  (EDVAR),A
          CALL SOREV
          CALL GETRRS
          JP  RRE1

50  OFFRCK: PUSH DE

          CALL SOREV
          CALL GETRRS
55  LD  DE,TIME1
          CALL CKZERO

```



POP DE  
RET

```
*****
5  *
  *   LOADS STATION DATA IN STANUB TO STADAT - REVERSE
  *
*****

10 GETRRS: PUSH AF
    PUSH BC
    PUSH HL
    PUSH DE

15    LD A,(STANUB)
    CP 16
    JP P,GSTR1
    LD HL,STARR
    DEC A
20    CALL GETPT
    LD DE,STADAT
    LD BC,31
    LDIR
    LD HL,STADAT
25    LD DE,STADAT1
    LD BC,31
    LDIR

    GSTR1: POP DE
30    POP HL
    POP BC
    POP AF
    RET

35 *****
  *
  *   SAVES STATION DATA IN STADAT TO STATION IN STANUB - REVERSE
  *
*****

40 SAVRRS: PUSH AF
    PUSH BC
    PUSH HL
    PUSH DE

45    LD A,(STANUB)
    CP 16
    JP P,SSTR1
    LD HL,STARR
50    DEC A
    CALL GETPT
    PUSH HL
    POP DE
    LD HL,STADAT
55    LD BC,31
    LDIR
    LD HL,STADAT
```

```

LD DE,STADAT1
LD BC,31
LDIR

5  SSTR1: POP DE
      POP HL
      POP BC
      POP AF
      RET

10 *****
   *
   *   EDIT SOLUTION LIST - REVERSE
   *
15 *****

REVSOL: LD A,83
      LD (SCRNUB),A
      CALL DISPSR

20      CALL ORSOL

REVSL1: LD A,(FLAG1)
      SET 7,A           ;REV VIDEO
25      LD (FLAG1),A

      LD A,(REVSTEP)
      CALL OUTSOL

30 REVSL2: CALL GETKEY
      CP 030H           ;SOLUTION LIST
      JP Z,RSOLST
      CP KBF1           ;DONE
      JP Z,REVRUN
35      CP KBCL          ;CLEAR
      JP Z,CLRSOL
      CP KBDN           ;FIELD ADV
      JP Z,RSADV
      CP KBUP           ;FIELD REV
40      JP Z,RSREV
      CP 0DH            ;ENTER
      JP Z,REVSL4
      AND 07FH
      SUB 030H
45      CP 0AH
      JP M,ENRSOL
      JR REVSL2

50 RSOLST: CALL SOLIST
      JR REV SOL

REVSL4: CALL SAVRRS
      JR REVSL1

55 CLRSOL: LD A,0
      LD (SOLUT1),A
      LD HL,PERCT1

```

```

        LD      A,20H
        CALL    FILVAR
        CALL    SAVRRS
        JR      REVSLI
5
ENRSOL:  CALL    REVIN
        CALL    SAVRRS
        JR      REVSLI

10 *****
   *
   *   ADVANCE TO NEXT STATION - REVERSE
   *
*****

15 RSADV: LD      A,(FLAG1)
        RES      7,A           ;NO REV VIDEO
        LD      (FLAG1),A

20      LD      A,(REVSTEP)
        CALL     OUTSOL

        INC     A
        CP      16
25      JR      NZ,RSAL
        LD      A,1
RSAL:   LD      (REVSTEP),A
        CALL    SOREV
        CALL    GETRRS
30      JP      REVSLI

*****
   *
   *   REVIEW PREVIOUS STATION - REVERSE
   *
35 *****

RSREV: LD      A,(FLAG1)
        RES      7,A           ;NO REV VIDEO
40      LD      (FLAG1),A

        LD      A,(REVSTEP)
        CALL     OUTSOL

45      DEC     A
        JR      NZ,RSR1
        LD      A,15
RSR1:   LD      (REVSTEP),A
        CALL    SOREV
50      CALL    GETRRS
        JP      REVSLI

*****
   *
55  *   OUTPUT 15 STATIONS WITH SOLUTIONS. - REVERSE
   *
*****

```

```
      ORSOL: PUSH HL
           PUSH BC

5      LD  A,1
      LD  (REVSTEP),A
      CALL SOREV
      CALL GETRRS

10     LD  A,(FLAG1)
      RES 7,A           ;NO REV VIDEO
      LD  (FLAG1),A

      RSOT1: LD  A,(REVSTEP)
15     CALL OUTSOL

           INC  A
           CP   16
           JR   Z,RSOT2
20     LD  (REVSTEP),A
      CALL SOREV
      CALL GETRRS
      JR   RSOT1

25     RSOT2 LD  A,1
           LD  (REVSTEP),A
           CALL SOREV
           CALL GETRRS

30     POP  BC
      POP  HL
      RET
```

```
*****
35    *
      *   STATION LIST - REVERSE
      *
*****
```

```
40    STARR: DW  RR1
           DW  RR2
           DW  RR3
           DW  RR4
           DW  RR5
45     DW  RR6
           DW  RR7
           DW  RR8
           DW  RR9
           DW  RR10
50     DW  RR11
           DW  RR12
           DW  RR13
           DW  RR14
           DW  RR15
```

```
55    *****
      *
```

\* STATION ORDER LIST - REVERSE  
\* GET STANUB FROM REVSTEP  
\*

\*\*\*\*\*

5

SOREV:PUSH AF  
PUSH BC  
PUSH HL

10

LD HL,RSTALST  
LD A,(REVSTEP)  
DEC A  
LD C,A  
LD B,0

15

ADD HL,BC  
LD A,(HL)  
LD (STANUB),A

20

POP HL  
POP BC  
POP AF  
RET

RSTALST DB 13  
25 DB 12  
DB 11  
DB 14  
DB 10  
DB 9  
30 DB 15  
DB 8  
DB 7  
DB 6  
DB 5  
35 DB 4  
DB 3  
DB 2  
DB 1

## CLAIMS

We claim:

1. An apparatus for automatically reprocessing a specimen from an infiltrating medium to an aqueous fluid comprising in combination:

a processing chamber for holding a specimen;

means for regulating flow of fluid to the processing chamber;

at least one container of a clearant agent, at least one container of a dehydrant agent and at least one container of an aqueous fluid, the containers of clearant, dehydrant and aqueous fluid being connected to the processing chamber via means for regulating flow of fluid to the processing chamber; and

a control device having a processor and a memory device, the processor controlling the means for regulating flow of fluid in order to automatically and sequentially connect the processing chamber with the container of clearant agent, the container of dehydrant agent and the container of aqueous solution in order to reprocess the specimen.

2. The apparatus of claim 1 wherein means for regulating flow of fluid includes a rotary valve and wherein the processor selects the containers of clearant, dehydrant or aqueous fluid by setting the rotary valve.

3. The apparatus of claim 1 further comprising:

at least one container of an infiltrating medium being connected to the processing chamber by a second valve and wherein the processor controls the second

valve.

4. The apparatus of claim 3 wherein the processor further controls the means for regulating flow of fluid and the second valve in order to automatically and sequentially, after the connection to the container of aqueous fluid, connect the processing chamber with the container of dehydrant agent, the container of clearant and the container of infiltrating medium in order to process the specimen.

5. The apparatus of claim 1 further comprising a container of purge dehydrant being connected to the processing chamber by the means for regulating flow of fluid, the processor controlling the means for regulating flow of fluid in order to automatically and sequentially connect the processing chamber with the container of clearant agent, the container of purge dehydrant, the container of dehydrant agent and the container of aqueous solution in order to reprocess the specimen.

6. The apparatus of claim 5 further comprising a container of purge clearant being connected to the processing chamber by the means for regulating flow of fluid, the processor controlling the means for regulating flow of fluid in order to automatically and sequentially connect the processing chamber with the container of purge clearant, the container of clearant agent, the container of purge dehydrant, the container of dehydrant agent and the container of aqueous solution in order to reprocess the specimen.

7. A computer readable storage medium containing a set of instructions

for a general purpose computer having a user interface comprising means for input, an output driver for connections to at least one valve, the valve being connected to at least one container of a clearant agent, at least one container of a dehydrant agent and at least one container of an aqueous fluid, the set of instructions comprising:

an initiating routine operatively associated with said user interface for permitting a user to initiate reprocessing via the means for input, said means for input being associated with a reprocessing program accessible to said computer;

a run routine for implementing said reprocessing program selected by the user, the reprocessing program controlling the output drive to the valve in order to automatically and sequentially connect the valve to the container of clearant agent, the container of dehydrant agent and the container of aqueous solution for reprocessing of a specimen.

8. The computer readable storage medium of claim 7 wherein the valve is connected to a container of purge dehydrant and wherein

the reprocessing program automatically and sequentially connects the valve to the container of purge dehydrant after connection of the valve to the container of clearant agent.

9. The computer readable storage medium of claim 8 wherein the valve is connected to a container of purge clearant and wherein

the reprocessing program automatically and sequentially connects the valve to the container of purge clearant before connection of the valve to the container of clearant agent.



10. Method for automatically reprocessing a specimen using a specimen reprocessing machine having processor for controlling the exposure of the specimen to a clearing agent, a dehydrating agent and an aqueous fluid, the method comprising the steps of:

providing the specimen which is infiltrated with an infiltrating medium;

indicating to the specimen reprocessing machine that the specimen is to be reprocessed;

exposing the specimen to a clearing agent via the processor to remove the infiltrating medium from the specimen; thereafter

exposing the specimen to a dehydrating agent via the processor to remove the clearing agent; and thereafter

exposing the specimen to an aqueous fluid via the processor to remove the dehydrating agent from the specimen.

11. The method of claim 10 wherein the processor further controls the exposure of the specimen to an infiltrating medium and further comprising the steps of:

exposing the specimen to a dehydrating agent via the processor after exposing the specimen to an aqueous fluid;

exposing the specimen to a clearing agent to remove the dehydrating agent;  
and

exposing the specimen to an infiltrating medium to replace the clearing agent.

12. The method of claim 10 wherein the clearing agent is xylene.
13. The method of claim 10 wherein the dehydrating agent is alcohol.
14. The method of claim 10 wherein the aqueous fluid is formalin.
15. A computer-readable storage medium containing a set of instructions for a general purpose computer, said set of instructions implementing the procedure shown in Figure 6a.
16. Method for reprocessing a specimen which is infiltrated with an infiltrating medium using a specimen reprocessing system having a processing chamber, the method comprising the steps of:
  - subjecting the specimen to at least one exposure to a clearant;
  - subjecting the tissue sample to at least one exposure to a purge dehydrant after the exposure to the clearant, the purge dehydrant being contaminated with clearant and being used to clean the processing chamber of the clearant;
  - subjecting the specimen to at least one exposure to a dehydrant after the exposure to the purge dehydrant; and
  - subjecting the specimen to an aqueous fluid.
17. The method of claim 16 further comprising the steps of:
  - subjecting the specimen to at least one exposure to paraffin; and thereafter
  - subjecting the specimen to at least one exposure to a purge clearant after the

exposure to the bath of paraffin, the purge clearant being contaminated with clearant and being used to clean the processing chamber of paraffin.

18. The method of claim 17 further comprising the steps of:

subjecting the specimen to the dehydrant after subjecting the specimen to an aqueous fluid;

subjecting the specimen to the clearant to remove the dehydrant; and

subjecting the specimen to an infiltrating medium to replace the clearant.

19. Method for reprocessing a specimen which is infiltrated with an infiltrating medium using a specimen reprocessing system having a processing chamber, the method comprising the steps of:

inputting a final step in reprocessing the specimen;

subjecting the specimen to at least one exposure to a clearant;

determining whether the step of subjecting the specimen to at least one exposure of clearant is the final step in reprocessing the specimen;

subjecting the specimen to at least one exposure to a dehydrant if the step of subjecting the specimen to at least one exposure of clearant is not the final step in reprocessing the specimen;

determining whether the step of subjecting the specimen to at least one exposure of dehydrant is the final step in reprocessing the specimen; and

subjecting the specimen to an aqueous fluid if the step of subjecting the specimen to at least one exposure of dehydrant is not the final step in reprocessing the specimen.

20. The method of claim 19 further comprising the step of

subjecting the specimen to an infiltrating medium to replace the clearant, after the step of determining whether the step of subjecting the specimen to at least one exposure of clearant is the final step, if the step of subjecting the specimen to at least one exposure of clearant is the final step in reprocessing the specimen.

21. The method of claim 19 further comprising the steps of:

subjecting the specimen to the clearant to remove the dehydrant, after the step of determining whether the step of subjecting the specimen to at least one exposure of dehydrant is the final step, if the step of subjecting the specimen to at least one exposure of dehydrant is the final step in reprocessing the specimen; and

subjecting the specimen to an infiltrating medium to replace the clearant.

22. The method of claim 19 further comprising the steps of:

subjecting the specimen to the dehydrant after subjecting the specimen to an aqueous fluid;

subjecting the specimen to the clearant to remove the dehydrant; and

subjecting the specimen to an infiltrating medium to replace the clearant.

1/10

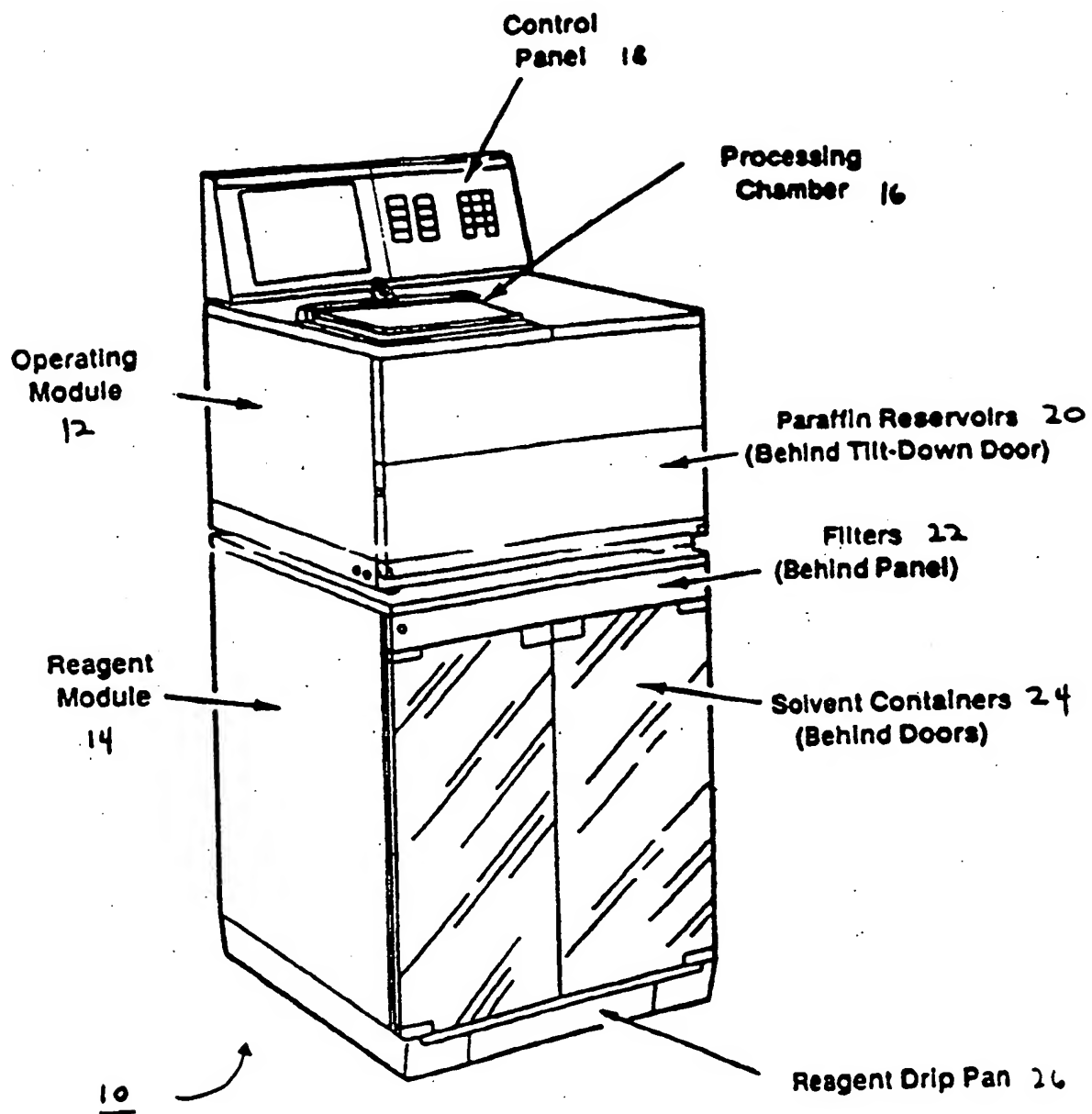


FIGURE 1

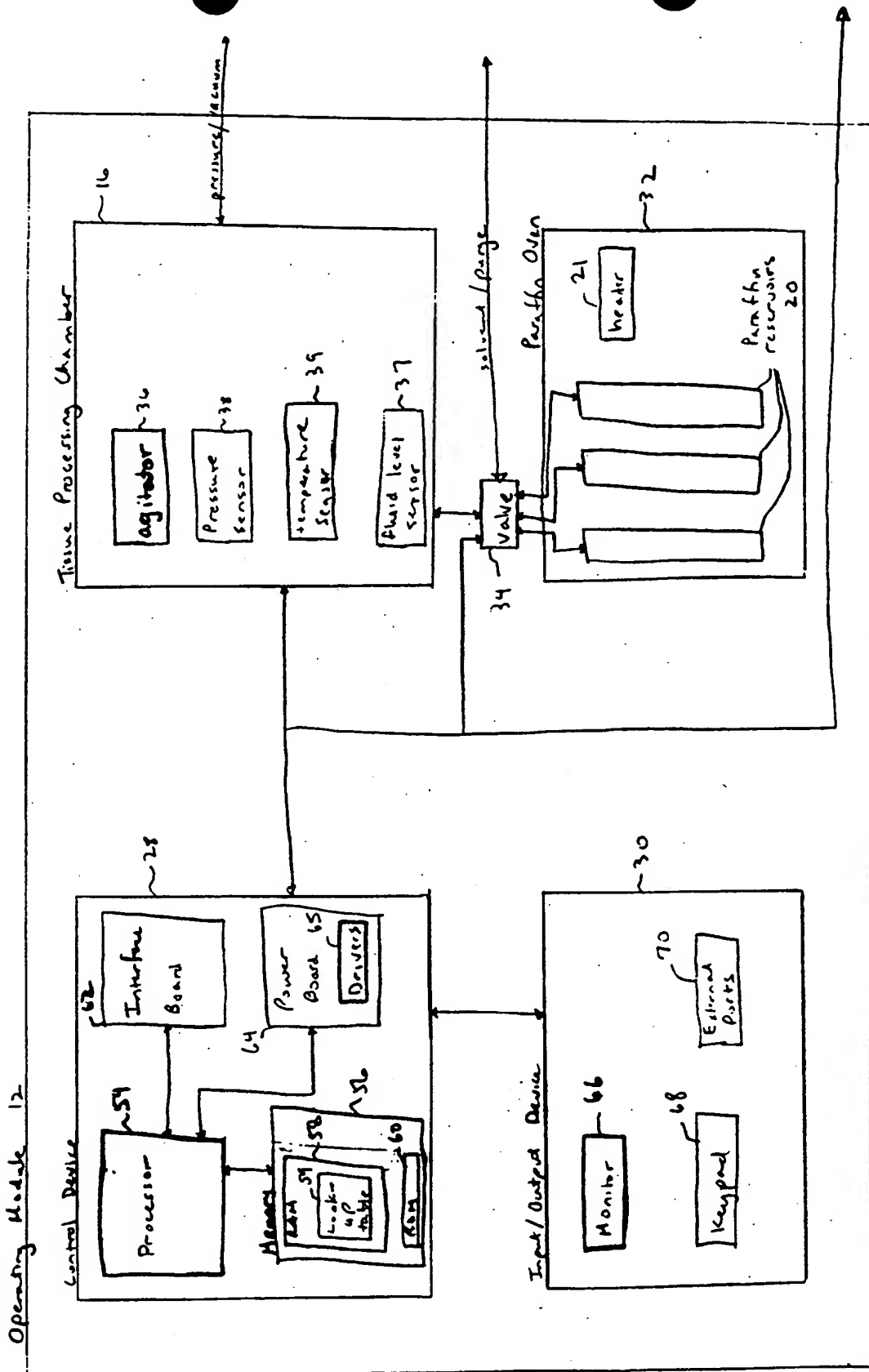


FIGURE 2a

3/10

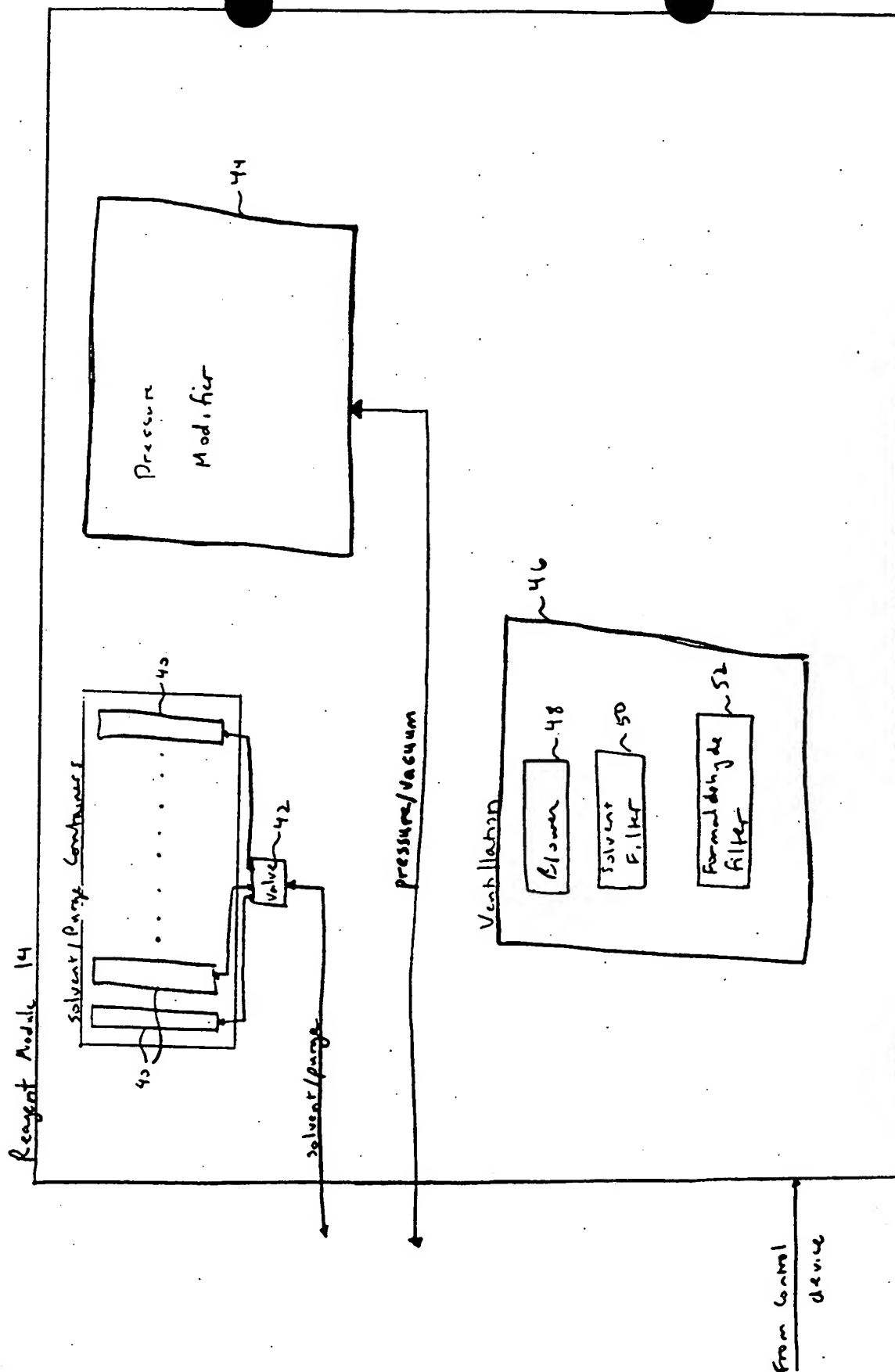


FIGURE 2B

4/10


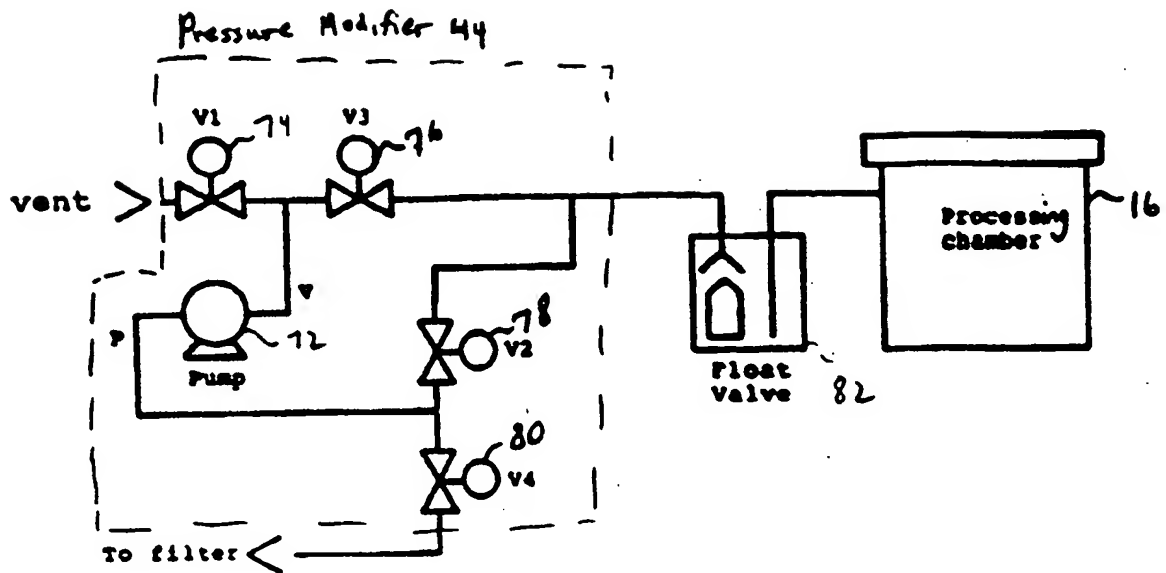


FIGURE 3



5/10

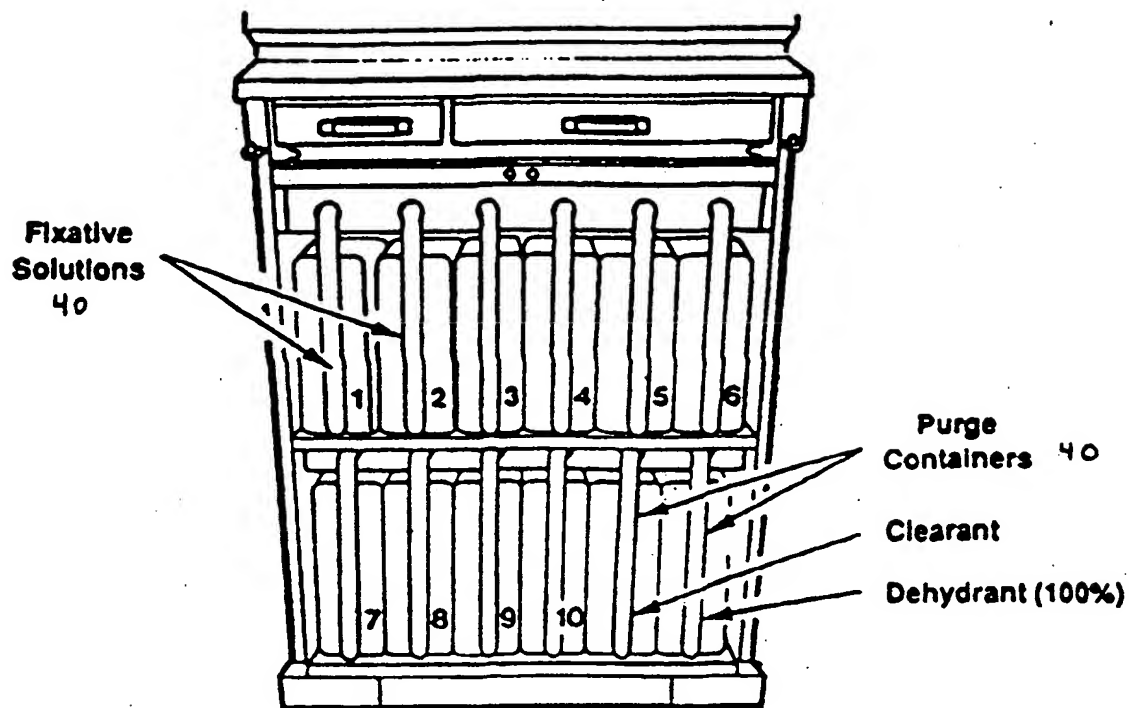

Solvent Container Positions in the Reagent Module

FIGURE 4

6/10

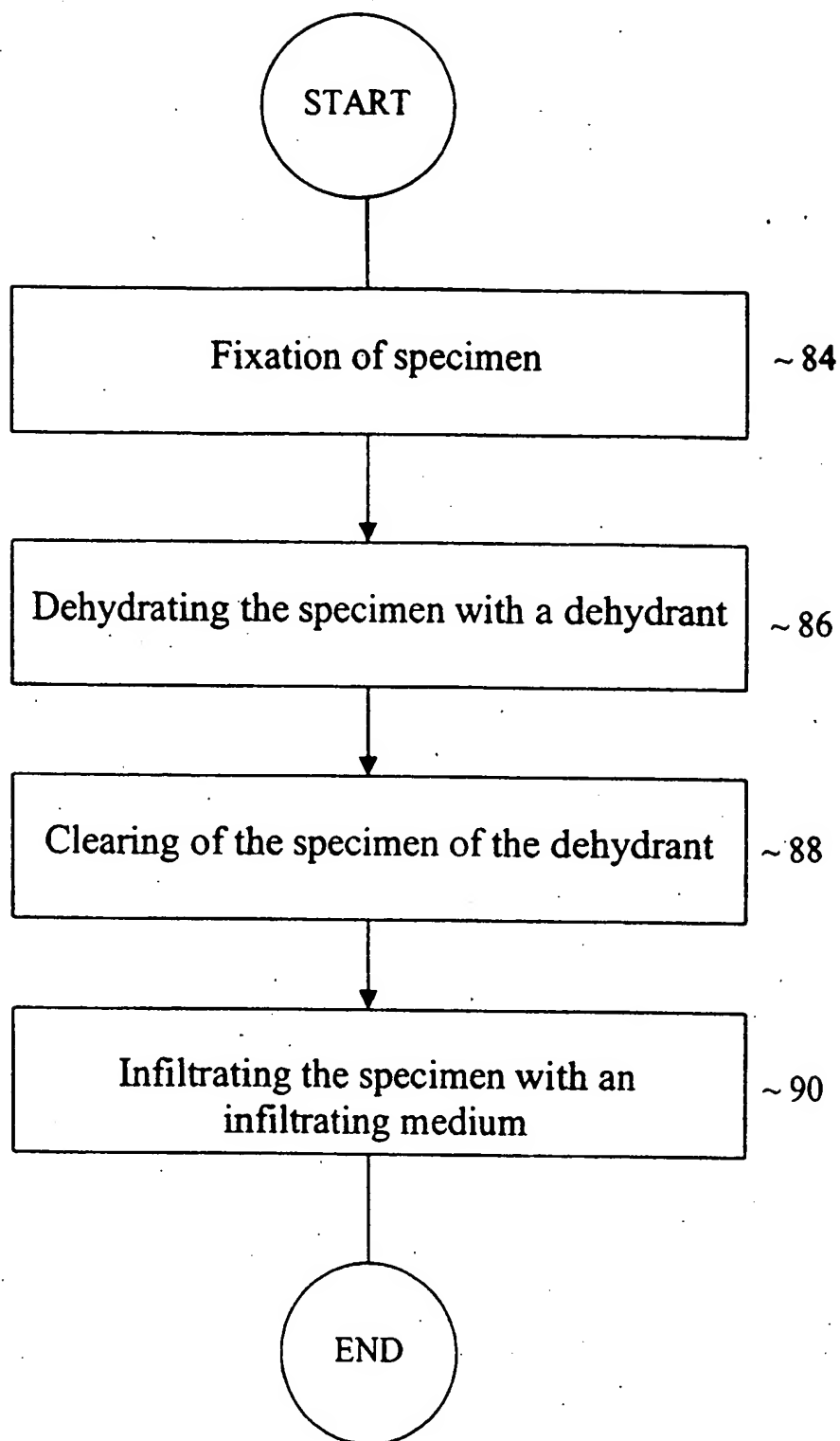


FIGURE 5

7/10


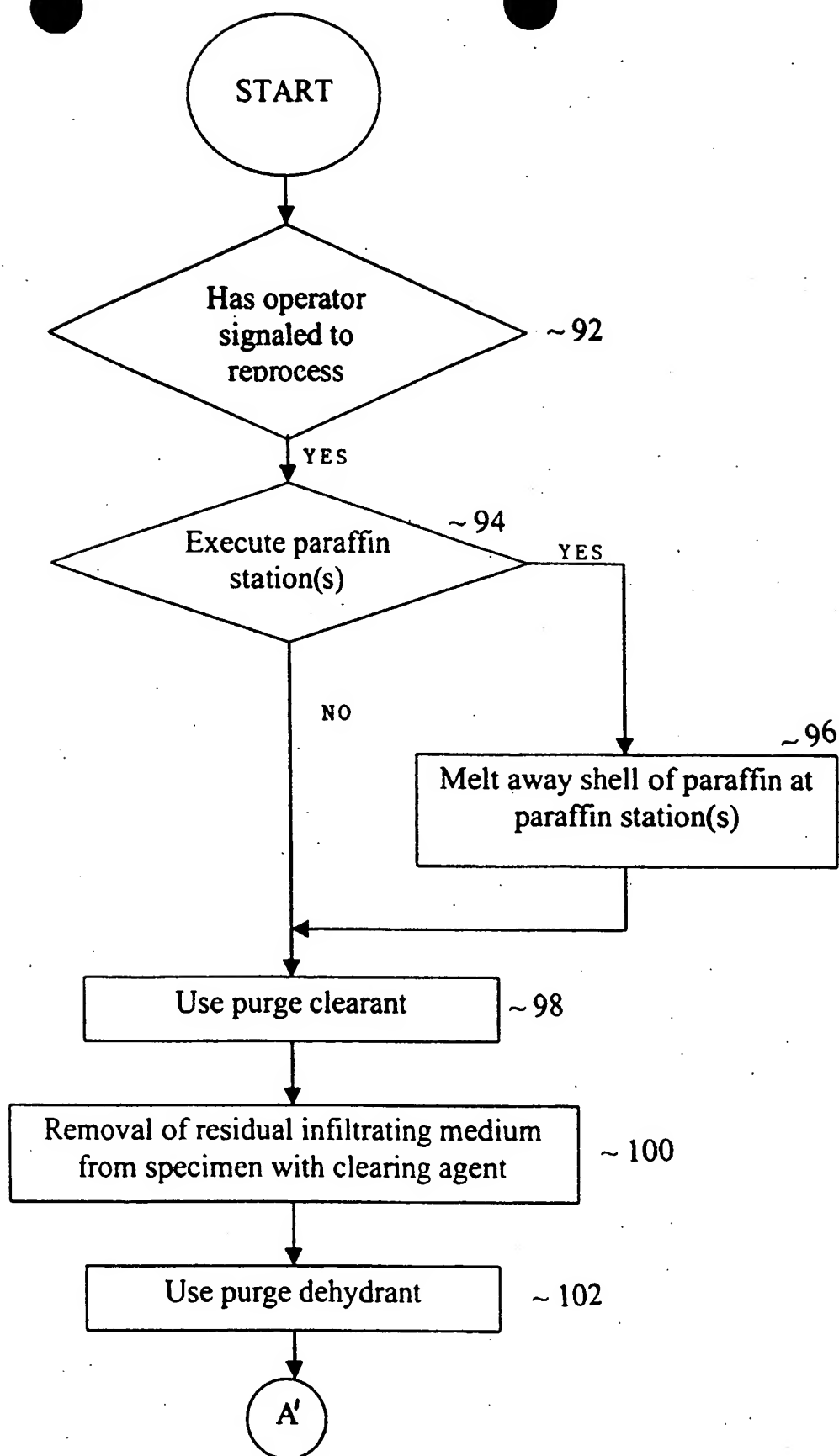


FIGURE 6a

8/10

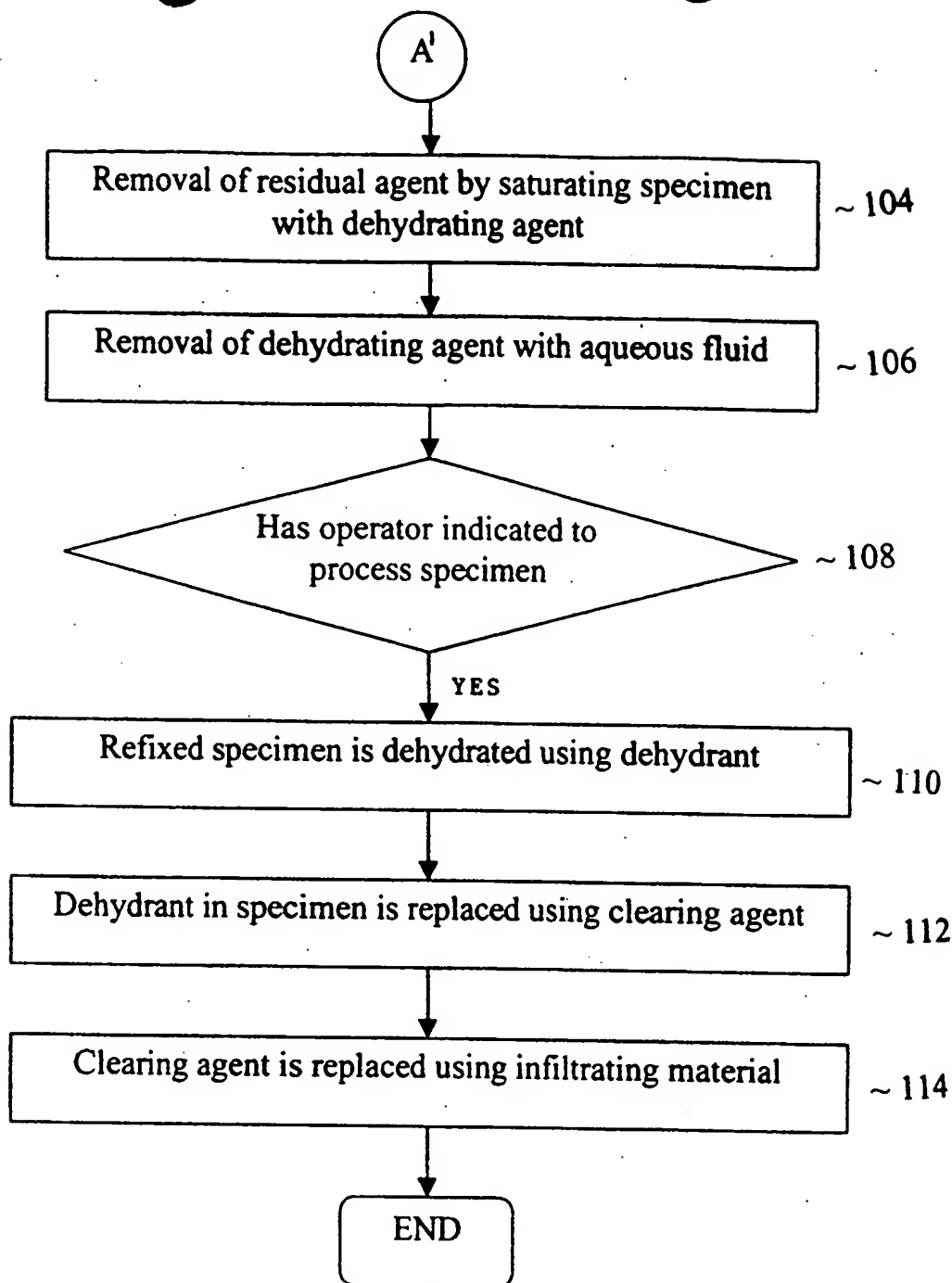

9/10

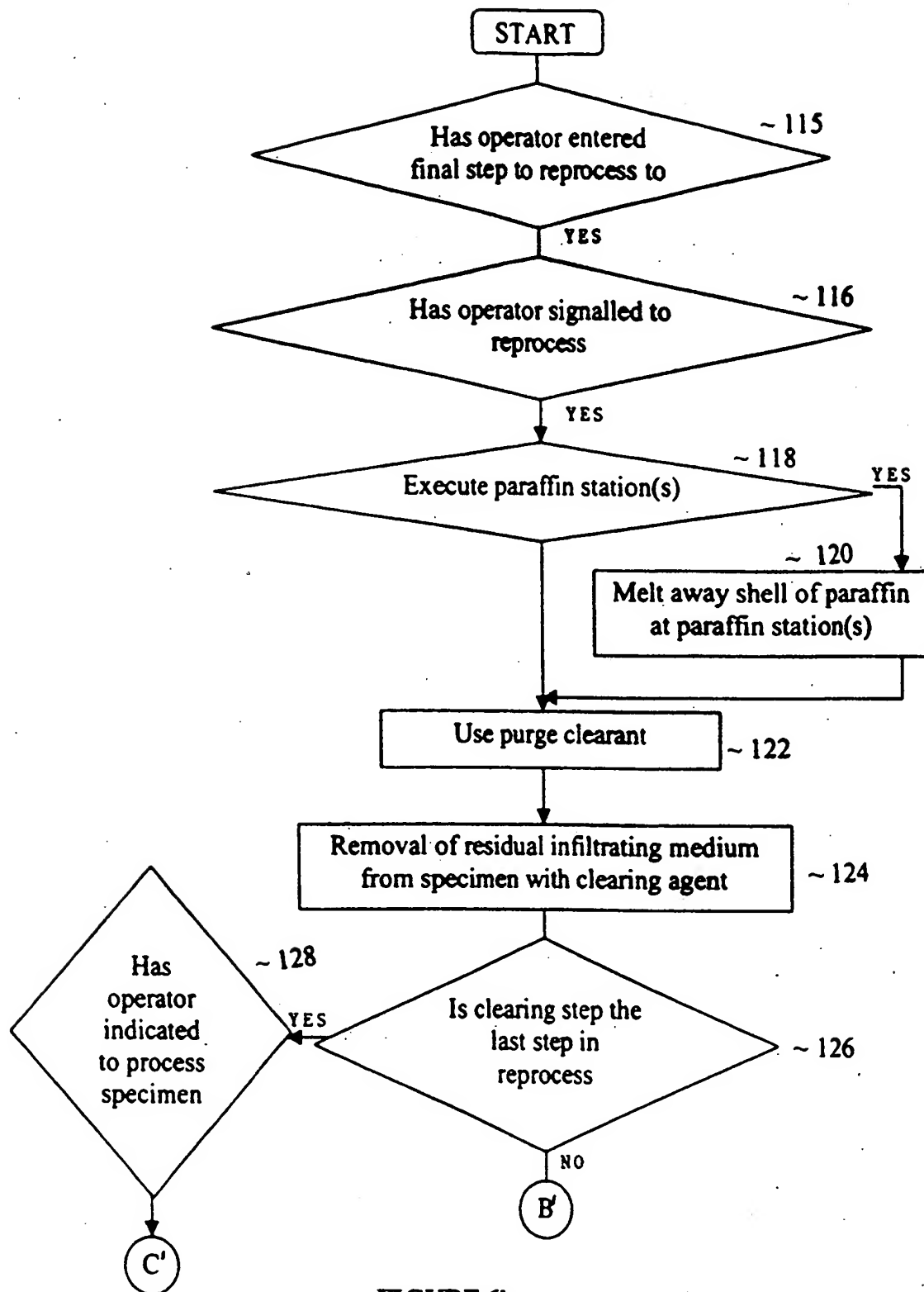


FIGURE 6b

10/10

